SKINERA™

Beautify your world

MICRODERMABRASION

- (EN) User Manual
- (DE) Bedienungsanleitung
- (FR) Mode d'emploi
- (NL) Gebruiksaanwijzing
- (ES) Manuale de usuario
- (IT) Manuale operativo

EN

DE

FR

NL

ES

IT

Table of contents

Get to know your Skinera Microdermabrasion	3
Contraindications & warnings	4
Do not use Skinera Microdermabrasion if:	7
What to expect with Skinera Microdermabrasion	7
Possible side effects.....	7
Before your first treatment	8
Get started!.....	10
Maintenance and storage	11
Spare parts and consumables	12
Troubleshooting.....	12
Customer service.....	13
Warranty	14
Disposal	17
Specifications	17
Signs and symbols.....	18

Get to know your Skinera Microdermabrasion

Skinera Microdermabrasion is an electrical device. As such, it should be used with special attention to safety. Before using the device for the first time, be sure to read all warnings, contraindications and safety precautions. It is important that you read and understand the instructions on how to use the device including after-use procedures, and follow them strictly. We recommend you refamiliarise yourself with the instructions before treatment.

Package contents

Your Skinera Microdermabrasion package includes these items:

- 1x Skinera Microdermabrasion device
- 1x Pack of filters
- 2x Treatment tips (Fine and Coarse)
- Spare rubber caps
- AC/DC adapter with cable
- User manual

Device description

Skinera Microdermabrasion is a dual-action therapy device designed for home use. It combines microdermabrasion exfoliation and vacuum stimulation of the skin.

Microdermabrasion exfoliation

- Gentle, non-abrasive exfoliation to allow for smoother, fresher skin.
- Helps improve skin tone, texture and elasticity.
- Helps reduce the appearance of age spots.
- Helps reduce the appearance of fine lines around the eyes.

Vacuum stimulation

- Helps stimulate for smoother, more youthful looking skin.
- Brings some impurities to the surface.
- Collects the dead skin cells for disposal.

The device has a treatment tip (#1), rubber cap (#2), filter (#3), control button (#4), device socket (#5), and an adapter with cable (#6). Refer to the diagram at the front of this manual to see the location of each feature on the device.

Intended use

Skinera Microdermabrasion is an over the counter device intended for the exfoliation of the skin's upper layer, for a fresher appearance. It is intended for the treatment of:

- Ageing skin of the face, neck and hands
- Fine facial lines
- Skin discolouration
- Dark spots and shallow scars
- Dull and sun-damaged skin

The device is intended for self-treatment at home. Do not use for any purpose other than as described in this manual.

Contraindications & warnings

Be sure to read this section carefully and take note of the warnings and contraindications.

This device must not be used by: children, people under the age of 16, those with limited physical, sensory or psychological capacities, or

people lacking experience or knowledge in how to use the device. Children should be supervised to ensure that they do not play with the device.

Do not use Skinera Microdermabrasion if you have:

- Rosacea.
- Poorly controlled endocrine disorders, such as diabetes and thyroid conditions.
- Abnormal skin conditions caused by diabetes, systemic or other metabolic diseases.
- Received radiation therapy or chemotherapy treatments within the past 3 months, or have areas of potential skin malignancies.
- Used physician-prescribed drugs on your skin within the past 6 months.
- Any other condition which, in the opinion of your physician or dermatologist, would make it unsafe for you to be treated.

Do not treat areas that have:

- Active eczema, psoriasis, lesions, open wounds or infections (cold sores), moles or warts, abnormal skin conditions caused by systemic or metabolic diseases (diabetes for example).
- A history of herpes outbreaks.
- Traumatized or injured skin.
- Had recent exposure to the sun without protection (wait 1 day); for tanned or sunburnt skin, wait 3 days.

Additional warnings

- Do not use on nipples or genitals.
- Recommended use is no more than twice a week.
- Use the low suction mode if you experience any discomfort.
- Stop use immediately if your skin is very red.
- Remove all jewellery before treatments.
- Protect your skin with SPF 30 or higher if you expose your skin to the sun after a treatment.
- Always ensure that your hands and skin are dry before operating the device or touching the adapter plug.
- Keep your device away from water or liquids! Do not use near bathtubs, showers, basins or other vessels containing water. Do not place or store the device where it can fall, be pushed or placed anywhere near water or any other liquid. There is a risk of electrocution. Unplug immediately and do not reuse.
- Never leave the device unattended when plugged into an outlet. Keep the AC/DC adapter accessible at all times, so it can be removed from the socket easily in case of emergency. Unplug the device immediately after use.
- Do not carry the device by the AC/DC adapter cable.
- If the device malfunctions, unplug immediately.
- Do not attempt to modify, open or repair your Skinera Microdermabrasion device. You may expose yourself to dangerous electrical components and may cause serious harm to yourself.

Do not use Skinera Microdermabrasion if:

- The adapter, attachments or accessories are not recommended by Invention Works B.V.
- The device has a damaged cable or plug. Keep the power cable away from heated surfaces.
- You see or smell smoke when it is in use.
- The device is not working properly, or it appears to be damaged.
- The device casing is cracked or coming apart.
- The device has fallen into water.

In such cases, contact our Customer Service.

What to expect with Skinera Microdermabrasion

During a treatment session it is normal to experience or feel:

- **A humming sound from the device:** This is the suction pump working.
- **Some slight redness of the skin:** This should subside within a few hours, or up to one day.

- Stop use immediately if your skin is very red.
- Contact your physician, if you experience any complications.

Possible side effects

If you follow the instructions fully and pay attention to all warnings and contraindications, complications and side effects are uncommon. All cosmetic procedures, including those designed for home use, involve some degree of risk. Therefore, it is important that you understand and accept the risks and complications that may occur.

Adverse reaction	Risk
<p>Skin redness</p> <p>You may experience dilation of the superficial blood capillaries. This may result in some redness for a few days. Additionally, some capillaries may present red spots that will also subside with time. In such cases, immediately.</p>	Minor
<p>Bruising</p> <p>You may experience slight bruising which should disappear after 2-3 days.</p>	Minor

Consult your physician if you experience any persistent symptoms.

Before your first treatment

To use your Skinera Microdermabrasion, follow these steps:

1. Remove the device and other components from the box.
2. Plug the adapter cable into the device socket and plug the adapter into an electrical outlet.
3. Review the recommendations on how to use Skinera Microdermabrasion on different treatment areas.

Treatment areas

We recommend that you start with the low suction mode and move to the high mode after a few uses. If you have any discomfort, work the treatment area in low suction mode.

Treatment area	Recommendations
Forehead, cheeks and chin	Move the device over your skin from the centre outwards and towards the hairline.

Eyes and nose	<p>Use short upward strokes for the area around your eyebrows. Avoid the eyelids.</p> <p>For the nose area, use straight downward strokes following the curve of the nose and the contour of the cheekbones.</p>
Neck	Use smooth, downward strokes and work the whole of the neck area.

To start with, you may work each area 2-4 times. When your skin is used to the treatment, you may increase to 4-6 passes.


Skinera Microdermabrasion's Fine and Coarse tips

The device has two treatment tips. We recommend that you start with the Fine tip.

- Use the Fine tip for general treatment or delicate finishing.
- Use the Coarse tip for heavy duty exfoliation.

To change tips:

1. Make sure the device is not connected to an electrical outlet.
2. Remove the current tip (#1) from the device and replace with the other tip. Make sure the tip is firmly in place before use.

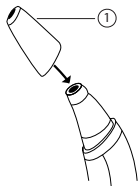Get started!

You're on your way to a fresher, smoother skin! Treating with Skinera Microdermabrasion is easy. Just follow these steps for each area you wish to treat. A total treatment should take around 5 minutes.

Recommended precautions

- Before you treat an area for the first time, do a spot test on a small area of skin. If you have no adverse reaction, go ahead and treat that area.
- Use the low suction mode for your first 2-3 treatments.
- We recommend you do not use the device more than twice a week.

- Stop immediately if you feel any discomfort.

Treatment steps

1. Make sure your skin is clean, dry and free from powder, cream or gel.
2. Make sure that the rubber cap (#2), filter (#3), and tip (#1) are firmly in place. For your first few treatments, be sure to use the Fine tip.
3. Press the control button (#4) to turn on the device. Your Skinera Microdermabrasion is ready for use in low suction mode.
4. Apply the tip firmly to your skin so that you feel the suction.
5. Work each area as recommended in section "Treatment areas".
6. To move from low to high suction mode, press the control button again.
7. To switch off the device after treatment, short press the control button once or twice, depending on which suction mode you are in.
8. Unplug the device and store Skinera Microdermabrasion in its original box. Please refer to section "Maintenance and storage".

After treatment

- Avoid direct sunlight immediately after treatment, and apply sun factor SPF 30 or higher.
- You may wish to apply a facial mask, or apply a day or night cream.

Maintenance and storage

It is important to clean your Skinera Microdermabrasion device and store it properly after each treatment.

- First, unplug the device.
- To clean the device, wipe with a clean, dry cloth. Use a cleaning fluid with no more than 70% alcohol, or another non-abrasive cleaning fluid. Never use acetone or other abrasive fluids, because they will damage your device.
- To clean the tip, first remove it from the device, wash with soap and water, then dry thoroughly.
- To clean the filter, first remove it from the device, brush away any debris using a soft brush, and place it back on the device.
- After cleaning, we recommend you store the device in its original box.
- With the exception of the tips, never immerse your Skinera Microdermabrasion in water.


Replacing the filter

The filter needs cleaning or replacing due to the build-up of discarded skin. We recommend that you do this every few treatments.

To replace the filter:

1. Remove the tip from the device. You will then see the rubber cap and the small black filter.

2. Remove the rubber cap (#2), and then remove the filter (#3).
3. To clean the filter, brush off any debris, then place it back on the device. To change a filter, remove the current filter and insert a new one.
4. Place the rubber cap back into position, then put the tip back on the device, making sure that it is firmly connected.

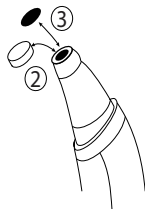Replacing the tip

We recommend that you replace the tip every six months, or earlier if the abrasive end of the tip is damaged.

Spare parts and consumables

You can order more filters and tips from our website: www.inno-webshop.com. For replacement rubber caps or AC/DC adapters, please contact our Customer Service.

Troubleshooting

Do not attempt to open or repair your Skinera Microdermabrasion. You may expose yourself to dangerous electrical components and cause serious harm to yourself. It will also void the warranty.

Problem	Checks
My Skinera Microdermabrasion does not start.	Make sure the adapter is properly connected to the device, that it is plugged into an electrical outlet and that you have pressed the control button.
I don't feel any suction when the device touches my skin.	Make sure the device is on, and check that the filter, rubber cap and tip are placed on the device correctly.

Only authorized personnel are permitted to perform repairs. For problems or queries, please contact our Customer Service.

- If the device fails, unplug immediately.

Customer service

See our website for more information: www.skinera.eu. This manual is also available as a PDF download from the website.

Contact Customer Service to report unexpected operations or events, or any other problem with your device. Our contact details are:

Country	Service number	Email
UK	0800-3586048	info@skinera.eu
Germany	04473-3270450	info@skinera.eu
France	0800-941298	info@skinera.eu

Netherlands	0900-0400492	info@skinera.eu
Spain	807001725	info@skinera.eu
Italy	-	info@skinera.eu

Warranty

What Does This Warranty Cover? Invention Works B.V., its applicable affiliates (hereinafter: Invention Works B.V.) or the authorized dealer from whom you purchased your Skinera™ product (hereinafter: the Device) shall repair or replace the Device, if any material defects caused by Invention Works B.V. are discovered in the Device during the warranty period, whether in the body of the Device or its operation.

Who is entitled to warranty by virtue of this warranty certificate?

Only you, the original buyer of the Device who purchased the Device from Invention Works B.V. or its authorized dealers, are entitled to receive the warranty specified in this warranty certificate. In order to be entitled to the warranty service, you need to provide a proof that you are the original buyer of the Device by presenting an original receipt for purchase of the Device. This warranty is not an international warranty. The warranty will only be honoured in the country in which the authorized dealer or the Invention Works B.V. entity from whom you purchased the Device, are located.

What period does this warranty certificate cover? The warranty covers repair requests received within one year (1) from the Device's original purchase date. Repair requests after one year (1) from the date of purchase are not covered by this certificate of warranty. Repair requests regarding a Device that was repaired or replaced must be submitted no later than: (a) one year (1) from the original purchase date of the Device or (b) ninety (90) days from the date that the repaired

or replaced Device was returned to you, the later to occur. If a sticker is attached to your Device stating a warranty period longer than that indicated above, the warranty period mentioned on the sticker will apply, and the Device will be repaired or replaced as part of the warranty until such period expires.

If you purchased the Device in an EU country, you have an additional year (according to the EU-Directive 1999/44) to claim for deficiencies or malfunctioning of your Device, which were present at the time of purchase, but were not apparent then, as determined by Invention Works B.V. or the authorized dealer. In order to make this extended warranty effective you must contact the authorized dealer from whom you purchased the Device or the local service center and notify them about the production fault immediately, and no later than 60 days after detecting the fault.

How to get service for the warranty? For further instructions how to obtain warranty service and for contact details of the service centres close to you, visit www.skinera.eu, or email info@skinera.eu.

When asking for customer services you may be asked to fill in a form for returning the product and you shall receive a reference number. In case of a defect that is not covered by the warranty for whatever reason, Invention Works B.V. can charge you for the shipping and repair costs if you agreed to repair the Device at your expense."

What is not covered by this certificate of warranty? This warranty does not cover any damage resulting from an accident, misuse (including using the Device other than as described in the use instructions accompanying the Device), defective use, lack of reasonable care, irregular and unreasonable use including without limitation commercial or professional use, adding parts not

supplied together with the Device, losing parts, using unauthorized spare parts or exposing the Device to any electrical current other than the indicated current. The warranty will be cancelled if the Device was treated by anyone not authorized by Invention Works B.V. and Invention Works B.V. will not be liable for any damage resulting from repairs performed by anyone not authorized by Invention Works B.V. The warranty does not include replacing used non-durable parts

WARRANTY DISCLAIMER THE EXPRESS. WARRANTY SET FORTH IN THIS DOCUMENT IS IN LIEU OF ALL OTHER WARRANTIES INCLUDING THE IMPLIED WARRANTIES OF MERCHANTABILITY, AND FITNESS FOR A PARTICULAR PURPOSE. THIS WARRANTY IS THE SOLE AND EXCLUSIVE REMEDY FOR ANY DEFECT OF THE DEVICE.

LIMITATIONS. IN NO EVENT SHALL INVENTION WORKS B.V., ITS AFFILIATES OR AUTHORIZED DEALERS, BE LIABLE FOR ANY INDIRECT, SPECIAL, PUNITIVE, EXEMPLARY, INCIDENTAL OR CONSEQUENTIAL DAMAGES (INCLUDING, WITHOUT LIMITATION, DAMAGES FOR LOSS OF USE OR CHARGES FOR SUBSTITUTE EQUIPMENT), REGARDLESS OF THE LEGAL THEORY ON WHICH THE CLAIM IS BASED EVEN IF THEY HAVE BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. IN THE EVENT THAT THE REMEDIES SET FORTH IN THIS WARRANTY ARE HELD TO FAIL OF THEIR ESSENTIAL PURPOSE OR IF INVENTION WORKS B.V., ITS AFFILIATES OR AUTHORIZED DEALERS ARE HELD TO BE LIABLE FOR DAMAGES ON THE BASIS OF ANY LEGAL THEORY, IN NO EVENT SHALL INVENTION WORKS B.V., ITS AFFILIATES OR AUTHORIZED DEALERS, LIABILITY FOR DAMAGES EXCEED THE WHOLESALE PRICE OF A SINGLE UNIT OF THE DEVICE AS OF THE DATE THAT THE CAUSE OF ACTION ACCRUES. YOU ASSUME ALL RISK AND LIABILITY FOR LOSS, DAMAGE OR INJURY TO YOU

AND YOUR PROPERTY AND TO OTHERS AND THEIR PROPERTY ARISING OUT OF THE USE OR MISUSE OF, OR INABILITY TO USE, THE DEVICE EXCEPT TO THE EXTENT CAUSED DIRECTLY BY THE NEGLIGENCE OF INVENTION WORKS B.V.

Disposal

Help protect the environment. Do not dispose with household waste and be sure to follow your country's laws for the disposal of electrical items.

Specifications

Model no.	IN4501
Technology	Air vacuum filling tip
Input power	12 V DC, 0.5 A
System weight	110 gm
Transport & storage between uses and storage condition	Temperature: -40°C to +70°C Relative humidity: 10% to 90% Atmospheric pressure: 500 to 1060 hPa
Operating conditions	Temperature: 10°C to 35°C Relative humidity: 30% to 75% Atmospheric pressure: 700 to 1060 hPa
Adapter models and rating	SAW12-050-2100XX XX Refer to plug type Input: 100-240VAC, 50-60 Hz, 0.5A Output: 12V, 0.5A

Expected service life of the equipment	3 years
Compliance	LVD: 2006/95/EC and 2014/35/EU EMC: 2004/108/EC and 2014/30/EU Directive 2011/65/EC RoHS REACH regulation EC No. 1907/2006

Signs and symbols

 Products marked with this symbol fulfil the requirements of the EU Directives.

Class III. The appliance is supplied from an extra-low voltage power source.

For indoor use only.

Do not use in wet environments.

Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE Directive.

Copyright © 2017 Invention Works B.V.
All rights reserved.

Invention Works B.V. reserves the right to modify its products or product specifications in order to improve their performance, reliability or manufacturability.

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Copyright 2017 © Invention Works B.V.